



## CourierPRO


### Passive Temperature Management Carrying Case

The Courier PRO is designed for healthcare providers and global immunization organizations requiring handheld portable solutions to transport vaccines, pharmaceuticals, healthcare specimens and other biomedical products.

The Courier PRO delivers reliable, long-lasting temperature control for up to 22 hours in a variety of temperature profiles (+2°C to +8°C, -10°C and -20°C) by using phase change coolant packs and vacuum insulated panel (VIP) technology. With the Courier PRO insulated carrying case, crucial temperature-sensitive healthcare products can be transported with confidence providing users with ultimate protection.



## Key Features

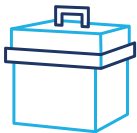

Lightweight, portable case



Superior VIP insulation, providing up to 10 times greater thermal efficacy




Available in three distinct temperature ranges



Phase change coolant packs provide long-lasting temperature control



Unique design eliminates need for ice or active power

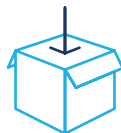

Simplistic conditioning and pack-out of coolant packs



Optimized design with 1 liter of payload space



Durable handle for comfortable carrying



**PAYLOAD (LITERS)**

**1**



**TEMPERATURE RANGE**

**+2°C to +8°C**

**-10°C**

**-20°C**



**DURATION (HOURS)**

**24-48**

MODEL	PAYLOAD	INTERNAL DIMENSIONS L x W x H	TARE WEIGHT	PACKOUT
Courier PRO	1 L 0.03 cu ft	236 x 84 x 50 mm 9.3 x 3.3 x 2 in	2.5 kg 5.5 lbs	Universal

MODEL	PRODUCT STORAGE TEMPERATURE	COOLANT PACK DIMENSIONS L x W x H	CONDITIONING METHOD	TOTAL WEIGHT
VXS-1 Coolant Packs	+2°C to +8°C	229 x 76 x 25 mm 9 x 3 x 1 in	Place in +3.5°C refrigerator for 24 hours	2.1 kg 4.6 lbs
VXS-2 Coolant Packs	at or below -10°C	229 x 76 x 25 mm 9 x 3 x 1 in	Place in conventional freezer for 48 hours	2.7 kg 6 lbs
VXS-4 Coolant Packs	at or below -20°C	229 x 76 x 25 mm 9 x 3 x 1 in	Place in conventional or lab/medical freezer at or below -28°C (-18.4°F) for 48 hours	2.1 kg 4.6 lbs

Contact us at [passivesolutions@csafeglobal.com](mailto:passivesolutions@csafeglobal.com)

CSafe has a worldwide network of comprehensive service centers, manufacturing facilities, sales offices and local hub sites to provide solutions where and when you need them. With representatives across the Americas, EMEA and APAC, local and regional support is always available.